


# Magnetic-Stirrer-MINI™

Popular small size MSM for CellVessel

Leaflet from CerCell who offer a stand-alone Magnetic-Stirrer-Mini (MSM) for adaptation to various small mixers and small SUBs. Suitable Working Volumes ranging up to 5 liter media.

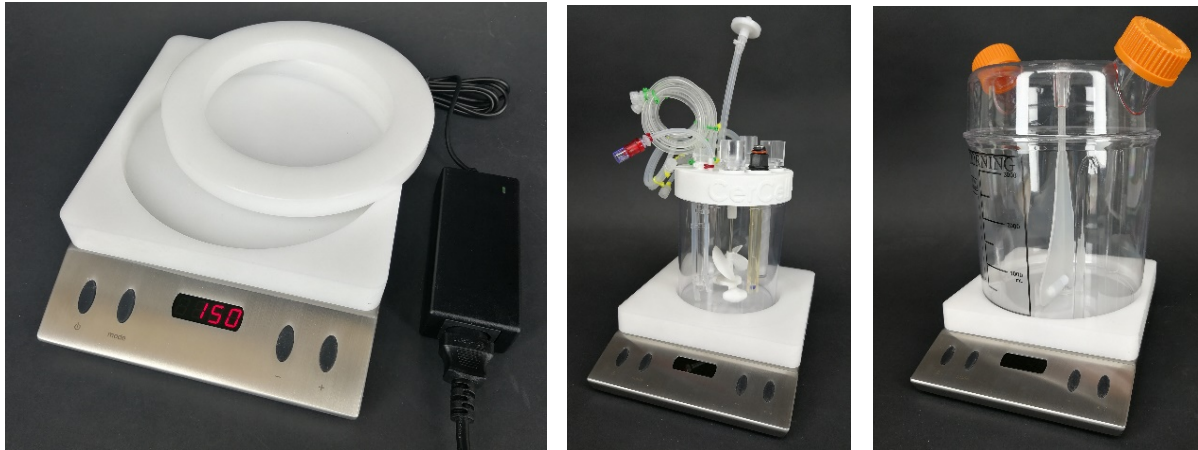

*All MSM facilitate a square block with ID 170 mm centring foot facilitating a range of inner-rings to support a variety of vessel diameters. Here shown mixers and SUBs in a stand-alone setup. Replaceable inner-rings ranging 82, 110, 120, 137 mm support various vessel diameters. Magnetic-Stirrer-Mini work equally well for spinner flasks for various applications. Able to handle the elevated temperature like 37°C in controlled environment.*

## Spec for MSM

- 100-2.000 RPM manually adjustable via 2 bottoms
- 15 watt max rotating / magnetic power transfer
- Material is stainless steel and fully encapsulated
- Max operating temperature +50°C
- Protection category is IP65
- Weight 4,5 kilo
- Supply voltage is 24 VDC
- Power supply CE marked and UL certified for 100-240 VAC with C13 power inlet
- Supplied with fixed ID 170 mm SSS and one extra centring ring either 82, 110, 120, 137 mm ID

September 2019, ver 1.1